

# INDIVIDUEL PRO

---

*by Individuel Genève*

## USER GUIDELINE

---

*Individuel*  
PRO  
GENÈVE

---


# INDEX

**Step 1** Unboxing

**Step 2** Assemble Technology Head & Capsule

**Step 3** Protocol: Capsule treatment

**Step 4:** After Usage Instructions

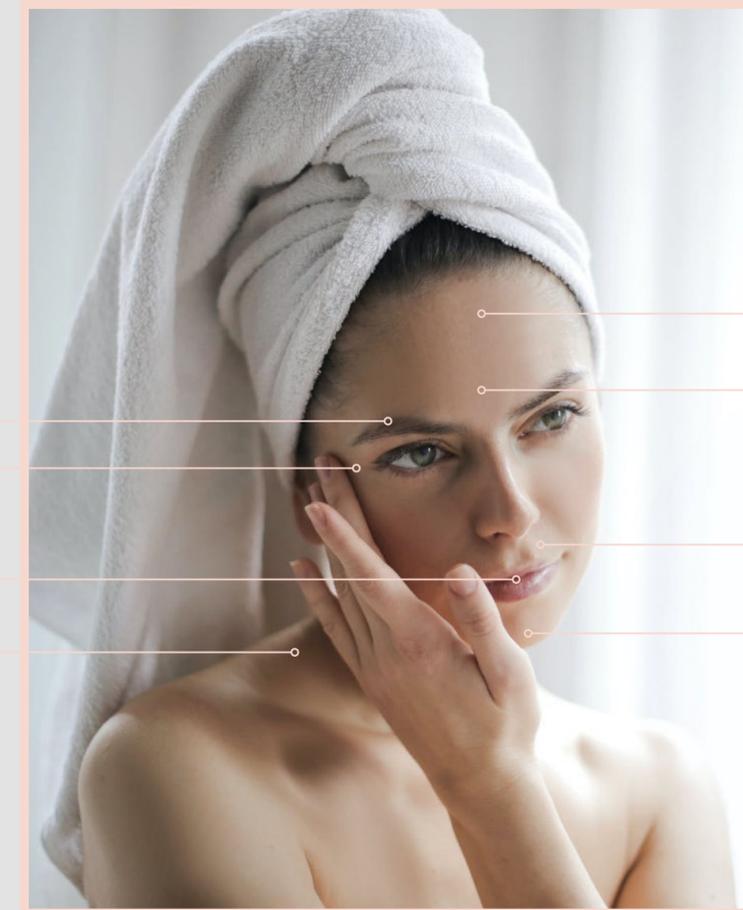

*Brow lift & asymmetric  
eyebrow shaping*

*Crow's feet  
(Peri-orbital lines)*

*Thinning lips*

*Neck folds  
)Platysma Muscle(*

*Forehead "worry" lines*

*Glabellar "Frown" lines*

*Wrinkles around lips &  
oral commissures*

*Puckered Chin*



# STEP 1

---

# UNBOXING

---

# UNBOXING

Inside:

- 1 Pink suitcase
- 2 Handle device
- 3 Charging base
- 4 Charger & Cord
- 5 2 Technology Heads
- 6 5 capsules
- 7 Folder with warranty and guideline



## STEP 2

---

### ASSEMBLE TECHNOLOGY HEAD & CAPSULE

---



#### IMPORTANT!

Be sure that your handpiece is fully charged prior to initial use.

Keep the handpiece at the charger whenever it is not in use.

1


# Connect Charger



2


Attach the technology head to the handpiece .

Gently **twist** the technology head clockwise until it clicks.

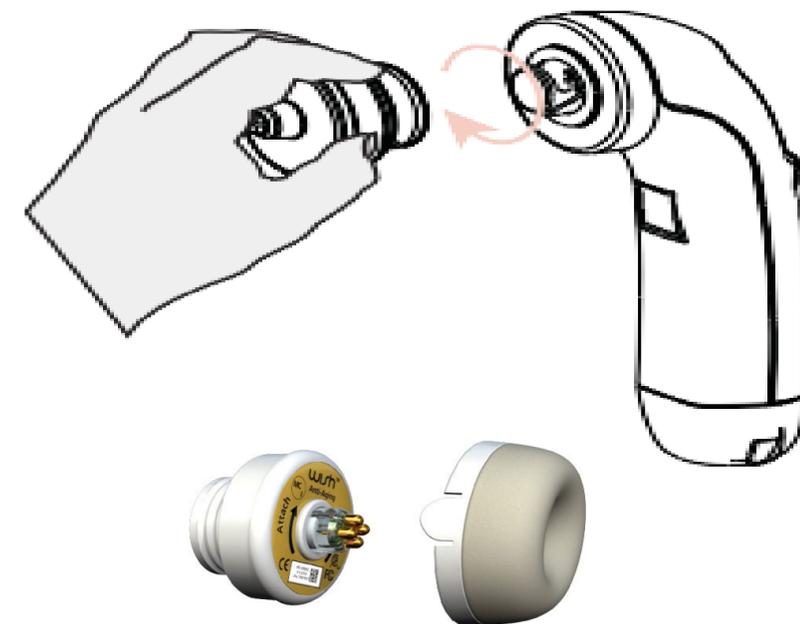

3

Hold the hand piece and **scan** the capsule against the technology head.

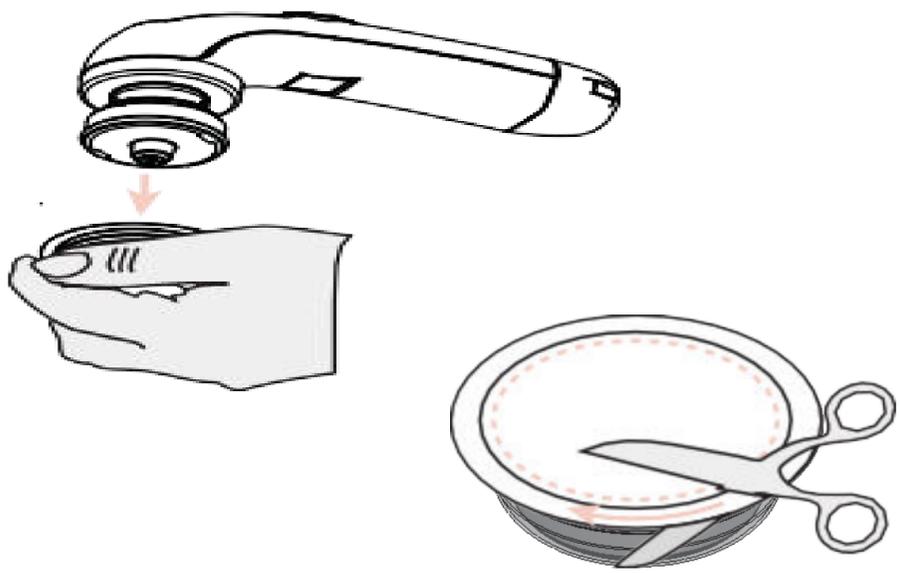

4



Open the appropriate capsule with a sharp tool .



5


Hold the capsules plastic and **attach** the capsule to the technology head.

Remove the plastic cup from capsule .  
(Keep the plastic cup until the end).

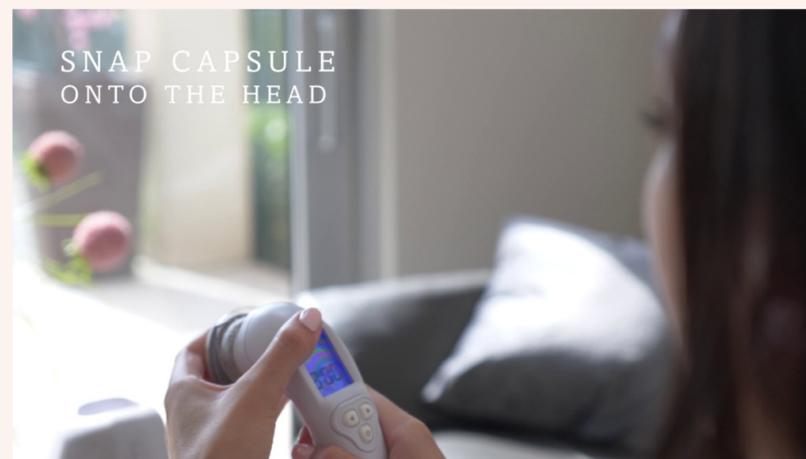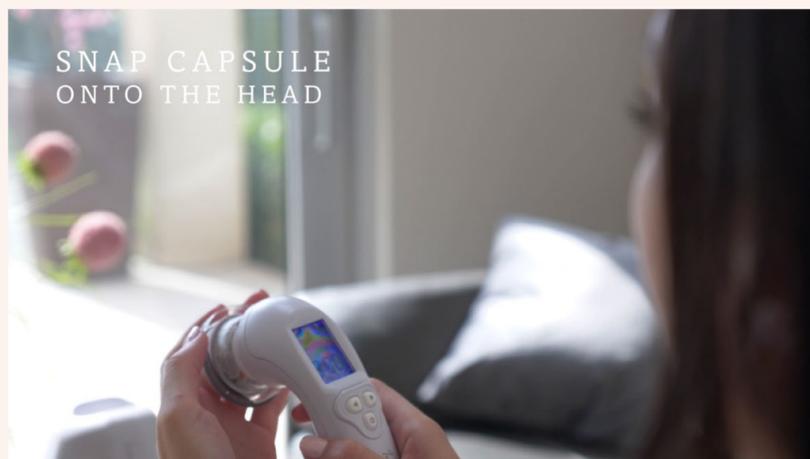


If capsule shows VOID after scanning DO NOT open it. Contact us at: [info@individuelcosmetics.com](mailto:info@individuelcosmetics.com) 0041 | 343 343 0800 | 78 86 810 22

5 Press the ON /OFF button to begin the procedure.

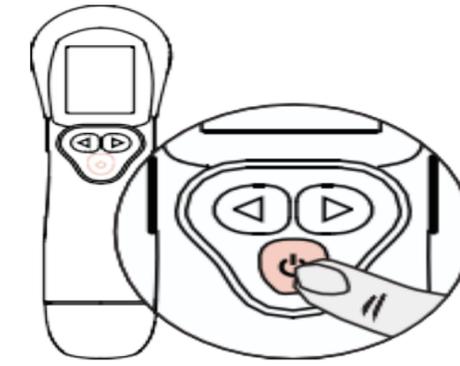

6 Adjust the intensity using the control buttons.  
Press left arrow twice to reduce intensity to half circle.

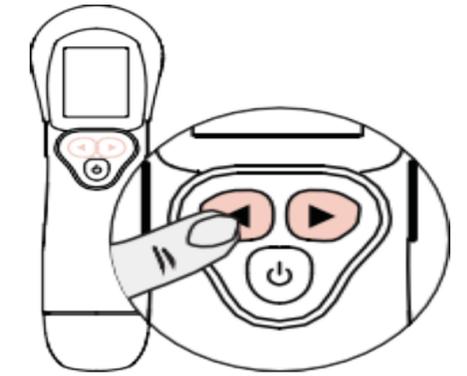


## STEP 3

---


### PROTOCOL: CAPSULE TREATMENT

---

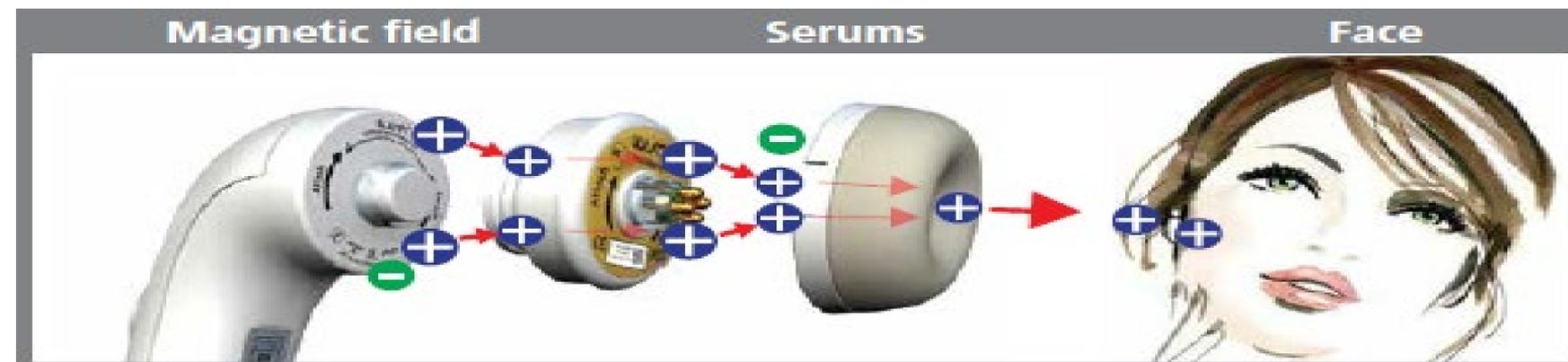

# PROTOCOL: CAPSULE TREATMENT

- 1 Clean your face with face cleanser
- 2 Place the head and the selected capsule.
- 3 Slide over the face & neck & décolleté pressing gently, with an upwards and outwards motion for 8 minutes until full absorbency.
- 4 When finished, do not wash face for at least 2 hours after the treatment.

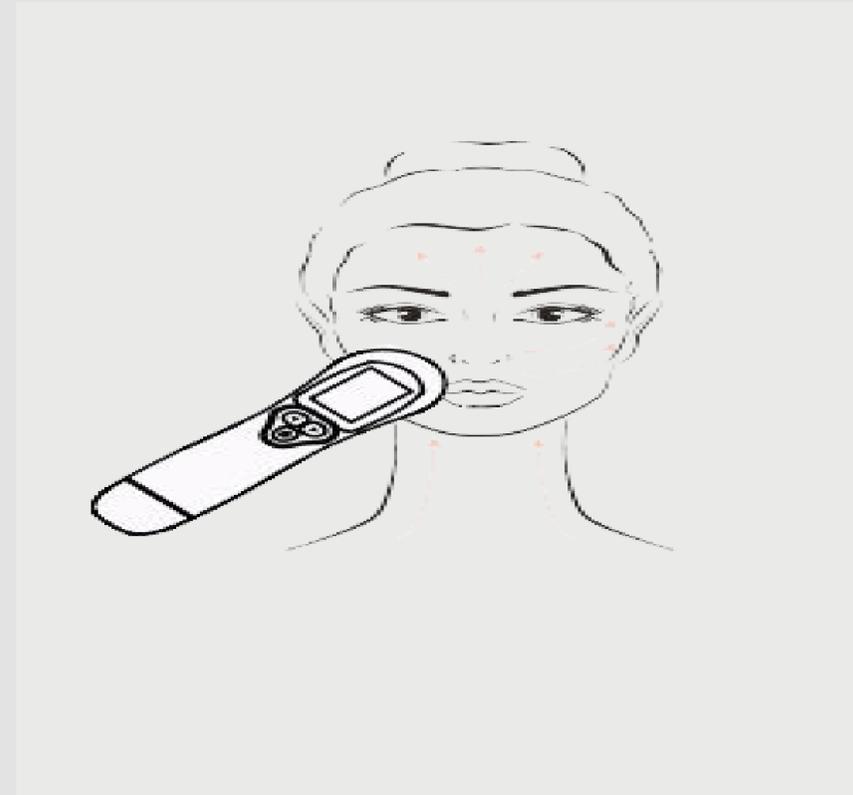



## STEP 4

---

### AFTER USAGE INSTRUCTIONS


---

#### IMPORTANT!


Do not grip the technology head by the rubber sleeve.

Be sure to grip the head on the hard plastic casing.

1 Shut off the device by pressing the ON/ OFF button.

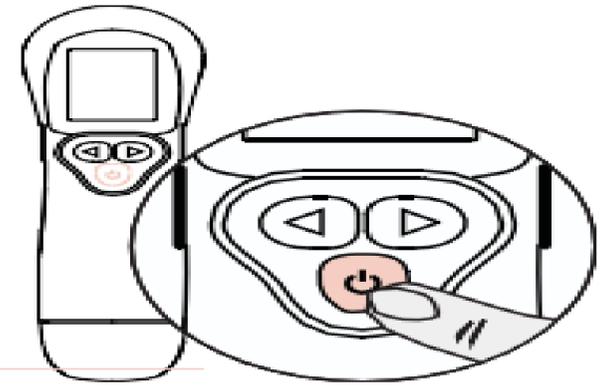


2 Firmly hold the handpiece in your right hand and remove the technology head with the capsule from hand piece by twisting it counter clock wise onto the plastic cup.

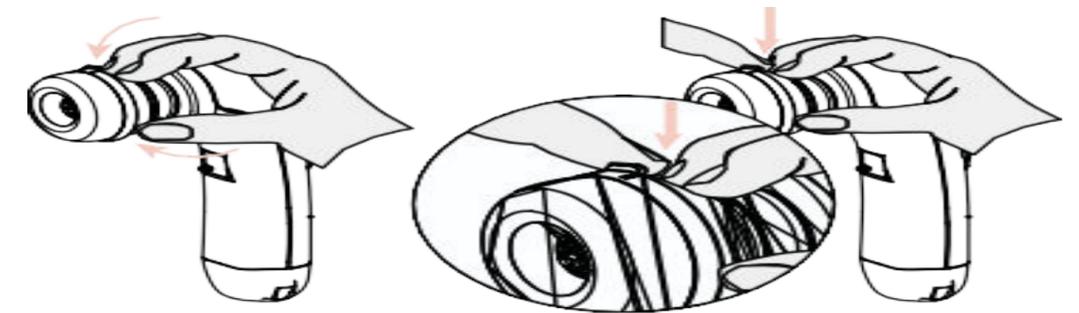


Insert the handpiece into the charger when not in use.

3

Detach the used capsule from the technology head and dispose it.

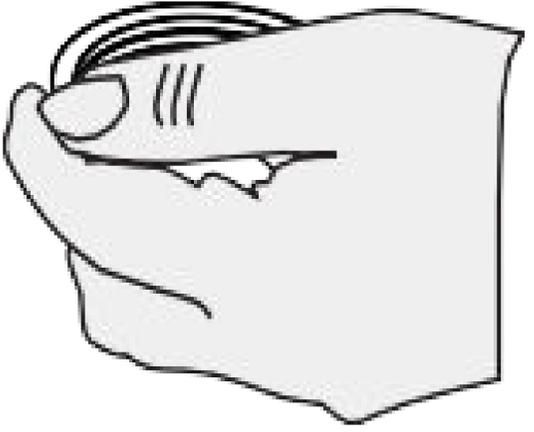

The capsules are designed for one- time use only.

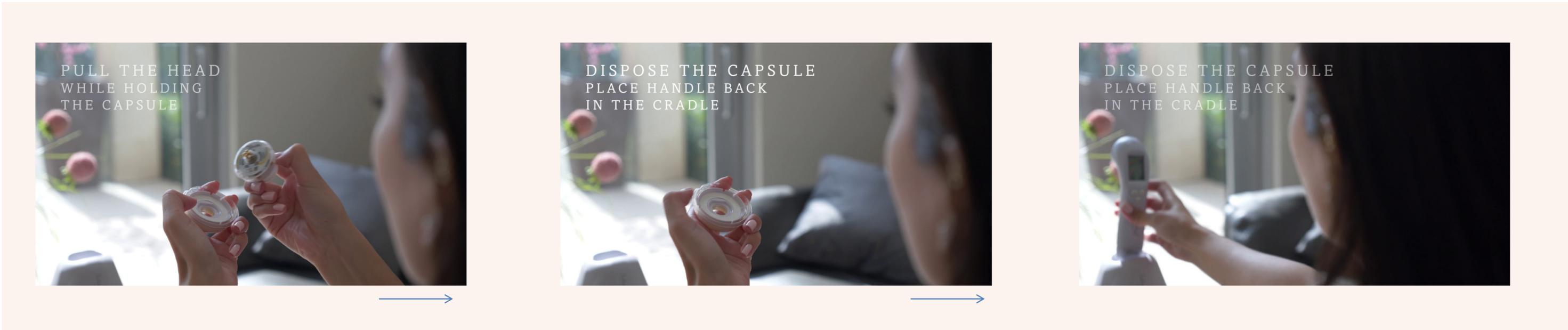

Insert the handpiece into the charger when not in use.

---

# INDIVIDUEL PRO

---

*Magnetic Field*    *Serum*    *Face*  
↓                    ↓                    ↓

